



YOUR GUIDE TO MULTI-ELECTRODE RENAL DENERVATION

The 8 F EnligHTN guiding catheter is designed to facilitate smooth device delivery and placement of the EnligHTN renal artery ablation catheter.

TRI-LAYER DESIGN FOR TORQUE, SUPPORT AND SMOOTH DELIVERY


- Multi-durometer polymer jacket
- Stainless steel braid
- Lubricious Polytetrafluoroethylene (PTFE) liner

HIGH VISIBILITY RDC-1 TIP

- RDC-1 designed for optimal renal positioning
- Radiopaque marker band 2 mm from the tip

HEMOSTATIC HUB AND SNAP-FIT DILATOR

- Hemostatic hub with a reliable seal design and snap-fit dilator¹
- Taper dilator tip for smooth insertion through introducer sheath
- 3-way stopcock with extension tubing

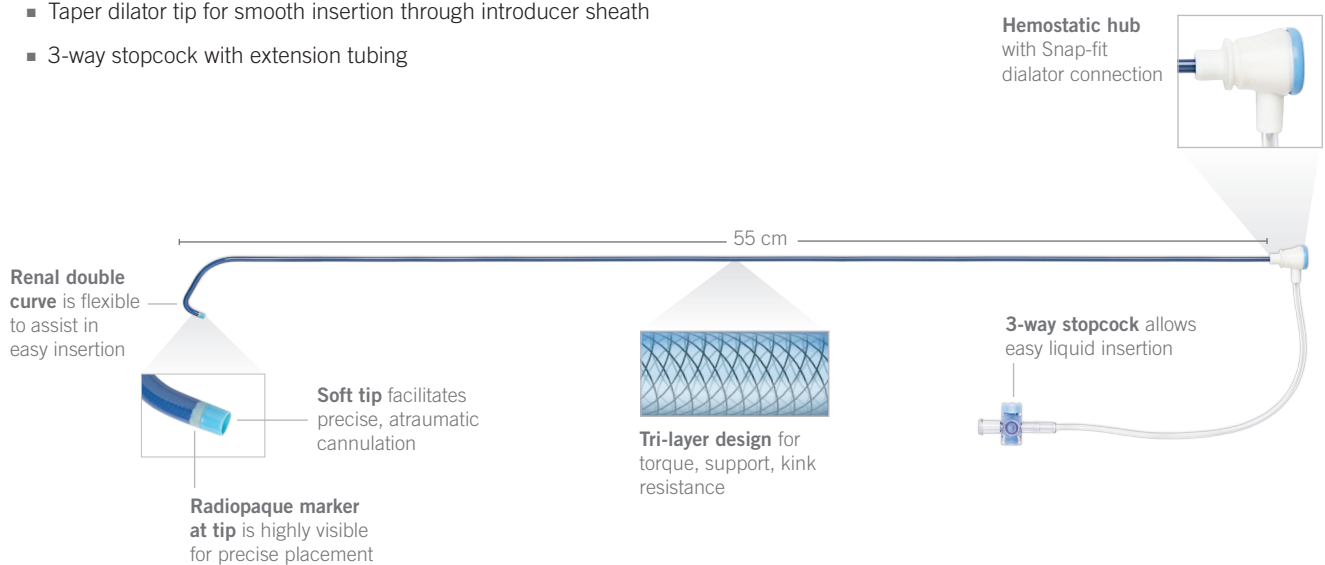

YOUR PARTNER IN RENAL DENERVATION TREATMENT

The EnligHTN renal denervation system offers advanced renal artery ablation treatment and is designed to deliver radiofrequency (RF) energy to the renal nerves to achieve targeted denervation. The core system components are the EnligHTN RF ablation generator, the EnligHTN renal artery ablation catheter and the EnligHTN guiding catheter.

Ordering Information

| EnligHTN Renal Denervation System | |
|-----------------------------------|---|
| Reorder Number | Product Name |
| ENL-8F55 | EnligHTN Guiding Catheter |
| ENL-LG-01 | EnligHTN Renal Artery Ablation Catheter - large |
| ENL-SM-01 | EnligHTN Renal Artery Ablation Catheter - small |
| ENL-T115 | EnligHTN RF Ablation Generator |
| ENL-CC-01 | EnligHTN 1641 Connector Cable |
| PREP-01 | Patient Return Electrode |
| ENL-PC-01 | EnligHTN Power Cable - replacement |



EnligHTN Guiding Catheter



EnligHTN RF Ablation Generator



EnligHTN Renal Artery Ablation Catheter

1. St. Jude Medical. Data on file.

ATRIAL FIBRILLATION

CARDIAC RHYTHM MANAGEMENT

CARDIOVASCULAR

NEUROMODULATION

Global Headquarters

One St. Jude Medical Drive
St. Paul, Minnesota 55117
USA
+1 651 756 2000
+1 651 756 3301 Fax

Cardiovascular Division

177 East County Road B
St. Paul, Minnesota 55117
USA
+1 651 756 4470
+1 651 756 4466 Fax

SJM Coordination Center BVBA

The Corporate Village
Da Vincilaan 11 Box F1
1935 Zaventem, Belgium
+32 2 774 68 11
+32 2 772 83 84 Fax

St. Jude Medical (Hong Kong) Ltd.

Suite 1608, 16/F Exchange Tower
33 Wang Chiu Road
Kowloon Bay, Kowloon
Hong Kong SAR
+852 2996 7688
+852 2956 0622 Fax

St. Jude Medical Brasil Ltda.

Rua Itapeva, 538
5º ao 8º andar
01332-000 – São Paulo – SP
Brazil
+55 11 5080 5400
+55 11 5080 5423 Fax

St. Jude Medical Japan Co., Ltd.

Shiodome City Center 15F
1-5-2, Higashi-Shinbashi
Minato-ku
Tokyo 105-7115
Japan
+81 3 6255 6370
+81 3 6255 6371 Fax

SJMenligHTN.com
SJMprofessional.com



ST. JUDE MEDICAL™
MORE CONTROL. LESS RISK.

Product referenced is approved for CE Mark.

EnligHTN, ST. JUDE MEDICAL, the nine-squares symbol and MORE CONTROL. LESS RISK. are registered and unregistered trademarks and service marks of St. Jude Medical, Inc. and its related companies. ©2012 St. Jude Medical, Inc. All rights reserved.

Item GMVASC402EN IPN 2157-12